



Tigsun COVID-19 Antigen Rapid Test

(Swab Sample by Lateral Flow Method)



Advantages

- Inexpensive, cost-efficiency ■
- Simple to use, no instrument needed ■
- Non-invasive, fast, result in 15 minutes ■
- Accurate and reliable compared to molecular methods ■

Fast / Convenient / Accurate / Reliable




● Principle


Tigsun COVID-19 Antigen Rapid Test employs immunochromatography technology to detect the SARS-CoV-2 antigen in human respiratory secretion specimens. Offering lab-quality results at the point of care with reliable performance in 15 minutes.

Reference :


Purpose and Options for Testing for SARS-Cov2(the COVID-19 virus):Considerations for World Bank Task Teams Managing COVID-19 Fast Track Facility Operations

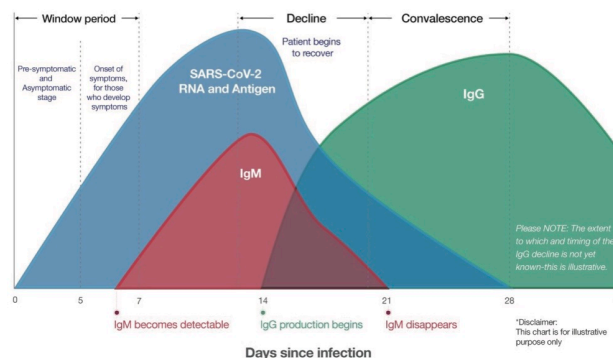

● Procedure

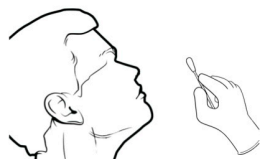


1 Respiratory secretion specimen collection

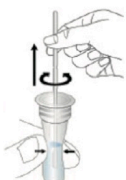

2 Add dilution buffer. Insert swab into reagent tube. Roll while pressing swab head



3 Close the tube with cap tightly

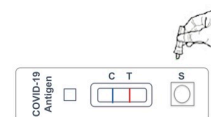

4 Read result within 15-20 mins after applying 3 drops to test card

● Test Result


Positive
(+)



Negative
(-)



Invalid

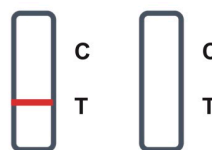

Product Name	Packing	Sample Type	Storage	Shelf Life
Tigsun COVID-19 Antigen Rapid Test	25T	Respiratory secretion specimen	2 ~ 30°C	24 months



Beijing Tigsun Diagnostics Co.,Ltd.

Add:No. 16 Guba Road, District #1, Chengguan Street, Fangshan District, Beijing, 102400, P.R.China.

Tel : +86-10-89350910 400-0666-093 (Toll Free)

Fax: +86-10-89350980

E-mail: world.sales@tigsun.com

Website: www.tigsun.com

